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Reagent kit for immuno-turbidimetric determination of Immunoglobulin G (IgG) in human serum.

Principle

Measurement of antigen/antibody reaction by the endpoint method.

Reference values

Serum: 800- 1600 mg/dl

It is recommended that each laboratory should assign its own normal range.

Reagents

1. Reagent 1 (R1)

buffer
sodium-azide (0,95 g/l)

2. Reagent 2 (R2)

antiserum
sodium-azide (0,95 g/l)

Precaution

Discard cloudy reagent. Avoid contamination by using clean laboratory material (pipettes, plastic vials for analyzers,...).

The reagents contain sodium azide. To avoid the possible build-up of azide compounds, flush waste-pipes with water after the disposal of undiluted reagent.

Sample

Use fresh serum. The serum can be stored at 2-8 °C for 48 hours. If stored for a longer period, the sample should be frozen.

Stability

without opening: till the expiry date indicated on the label
after opening: 14 days
calibration frequency: 7 days
onboard stability: 7-14 days
Stability data are valid only when using new system bottle!

PROCEDURE

The reagents are ready for use.

Assay conditions

Wavelength: 340 nm
Temperature: 37 °C
Cuvette: 1 cm light path
Read against: sample blank
Method: endpoint (increasing)

Pipette into cuvette

	Blank	Standard	Sample
1. reagent (R1)	360 µl	360 µl	360 µl
Standard	3 µl or 3 µl	3 µl	
Sample			3 µl

Mix, wait 20 seconds then add:

2. reagent (R2)		90 µl	90 µl
-----------------	--	-------	-------

Mix and after a 4-minute incubation read the absorbance against the sample blank.

Calibration

5 level (at the least 3) standard set for calibration is necessary.

S1: Distilled water or physiological salt solution
S2-S6: 5 level standard set

New calibration is recommended:

- after reagent lot change,
- as required following quality control procedures.

Calculation

Based on the calibration curve.

Quality control

A quality control program is recommended for all clinical laboratories. The analysis of control material in both the normal and abnormal ranges with each assay is recommended for monitoring the performance of the procedure. Each laboratory should establish corrective measures to be taken if values fall outside the limits. Recommended controls: Diagnosticum DunaCont N (Cat. No.: Dcon-N) and DunaCont P (Cat. No.: Dcon-P)

PERFORMANCES DATA

The following data were obtained using the Olympus AU600 and Advia 1650 analysers (37°C).

Linearity

The test is linear up to 3000 mg/dl (30 g/l).

Sensitivity

It is recommended that each laboratory establishes its own range of sensitivity as this is limited by the sensitivity of the spectrophotometer used.

The estimated value of limit of detection (LOD) is 3,7 mg/dl.

Precision

Sample	Reproducibility		
	Average concentration (g/l)	SD	CV %
sample I	10,17	0,53	5,23
sample II	9,05	0,63	7,01

Sample	Repeatability		
	Average concentration (g/l)	SD	CV%
sample I.	18,74	0,27	1,45
sample II.	11,51	0,04	0,31
sample III.	6,40	0,05	0,77

Correlation

Comparative studies were done to compare our reagent with another commercial IgG reagent. The results from these studies are detailed below.

Correlation coefficient: $r=0,9947$

Linear regression: $y (g/l) = 0,995x + 0,535$

(x= other commercial reagent, y= own reagent).

Specificity

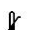
Bilirubin 20 mg/dl, triglycerides 1000 mg/dl, hemoglobin 1000 mg/dl don't interfere with the assay up to the given levels.

NOTE

Do not use reagents after the expiry date stated on each reagent container label. Do not use products, test solutions and reagents described above for any purpose other than described herein.

For *in vitro* diagnostic use only.

The following symbols can be used on the labels

 In vitro diagnostic device	 Batch code
 Manufacturer	 Catalogue number
 CE-marking	 This way up
 Temperature limitations	 Biological risk
 Use by (year/month)	

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